

# MYOPIA APP

## USER MANUAL

# Copyright Statement

© Copyright 2020 MDBootstrap.com

Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions. The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL

THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT, TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

## CONTENTS

01. INTRODUCTION .....	3
1.1 What is myopia app? .....	3
1.2 How to install myopia app.....	4
1.3 How we ensure the security of your data .....	5
02. PROGRAM workspace .....	5
2.1 Signing in to your account .....	5
2.2 App overview .....	6
03. How to .....	9
3.1 how to do the test.....	9
3.2 How to clarify the results.....	10

## 01. INTRODUCTION




### 1.1 WHAT IS MYOPIA APP?

Myopia app is a mobile application specially has been made up for myopia patients. This app allows you to measure the **visual acuity** of your vision by using your mobile phone or a tablet to get easy access to this anytime and anywhere.

#### Which devices does the mobile app support?

You can install Myopia app on any mobile devices that have Android operating systems.

### Responsive For All

	PC WEB
	IPAD/TAB
	MOBILE

#### How it works

Myopia (nearsightedness) is a refractive error. Refractive error is when the eye does not bend (refract) light properly. Light does not focus correctly so images are not clear. In myopia, close objects look clear but distant objects appear blurred. Myopia is a common condition that affects to whole world. It is an eye focusing disorder, not an eye disease. Now this disorder can easy identify using our web application. A user who is use this application he/she can save their cost and time. Further, this web application is provide medical and technical support for doctors and eye-specialist. Before use this web application we recommend some requirements should follow by the user.

#### Requirements

- Full HD (1920x1080) 15 inch monitor.
- 20ft or 10ft distance between user and the device.

- This testing should do in the normal lighted environment.

Device and the user's eyes should in the same level.

This is a web application which has include an electronic Snellen chart within it. Using this electronic Snellen chart, particular UI has been developed to read their characters to users. While checking visual acuity, mobile application always maintains the brightness of the display and distance between electronic Snellen charts and users. Further App showing characters of the Snellen chart on the display and the user has to pronounce them. The letter or number on the screen and user's answer is recorded inside the app. After testing the user's eyes, app provides a UI which shows the accuracy of the user's predictions. This measurement is used to measure the visual acuity of the peoples' vision in this app.

#### **Myopia app allows you:**

- To measure the eye sight level by yourself
- To compare the results with previous readings provided by the app
- To capture the users voice command accurately.

#### **How much does the mobile app cost?**

Currently, we only offer monthly subscription. You can upgrade or cancel your monthly account at any time with no further obligation. You can cancel your subscription anytime in your account. Once the subscription is cancelled, you will not be charged next month. You will continue to have access to your account until your current subscription expires. Unfortunately, we do not issue full or partial refunds for any reason.

## **1.2 HOW TO INSTALL MYOPIA APP**

Depending on your mobile device, go to Google Play and search for the Myopia app.

### **To find and install Myopia app for android:**

1. On your computer device, open **web browser**.
2. Tap the Search icon.
3. Enter Myopia app in the search field.
4. Select Myopia app in the search results to go to the app page.
5. Follow the standard installation procedure.

## **1.3 HOW WE ENSURE THE SECURITY OF YOUR DATA**

We work with top payment companies which guarantees your safety and security. All billing information is stored on our payment processing partner which has the most stringent level of certification available in the payments industry. Furthermore we store your data on our servers in encrypted form. Only you have access to your encrypted data.

## **02. PROGRAM WORKSPACE**

### **2.1 SIGNING IN TO YOUR ACCOUNT**

To use the web app, you need to sign in using the email address and password of your Myopia app account.

If you do not have a Myopia app account yet, click “Register” and follow the instructions.

Welcome to APP

Log in

E-mail

Password

[Forgot password?](#)

LOG IN

Not a member? [Register](#)




The Myopia app interface consists of three tabs-Home, Profile, and Log Out

Myopia APP   Home   Profile   Log Out

## Home Tab

By scrolling down in the home tab mainly, you can find “start” button which allows you to start the eye test and two other buttons which allows to have quick access to two other tabs respectively. Also you can find some important information about this app by scrolling down in this home tab.

Click the “Start” Button to start your eye test.

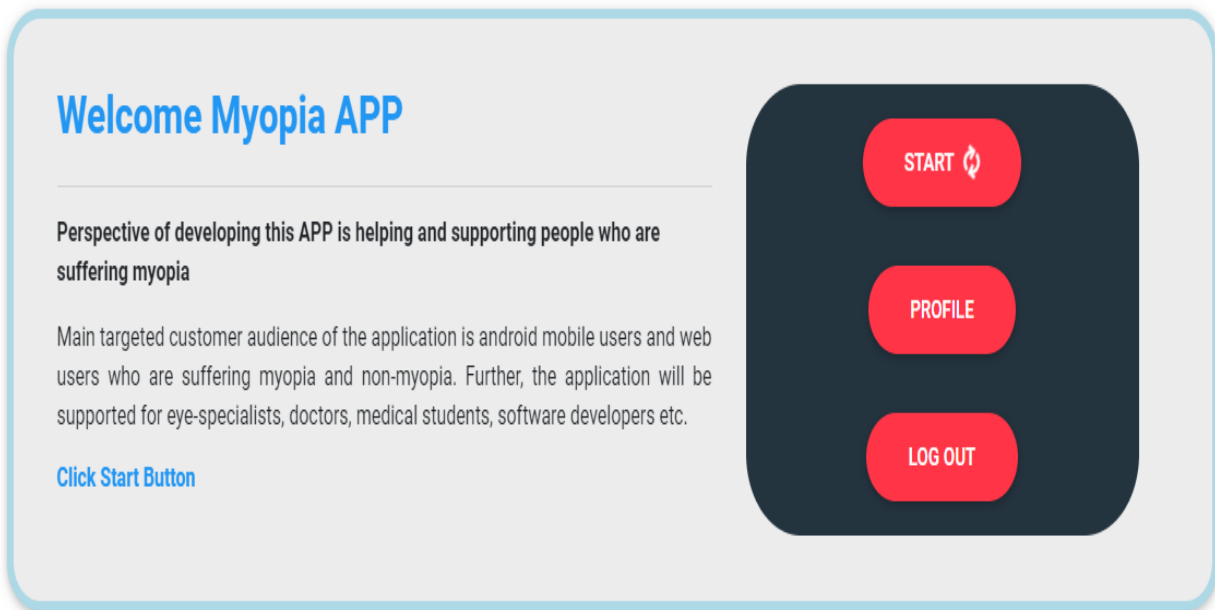

## Profile Tab

Click the profile tab or “Profile” button in the home tab to view your profile.

### Profile

First Name

Last Name

E-mail

Age

Gender

SHOW PAST REPORTS

By clicking the “Show Past Reports” button in the profile tab you can view your past reports as follow,

Myopia APP

Home

Profile

Log Out

{username}'s Reports

Test ID	Date & Time	Right Eye %	Left Eye %
1	20/11/2020	90%	50%
2	20/11/2020	30%	80%
3	20/11/2020	40%	10%

f in @

© 2020 Copyright: MDBootstrap.com

## Log

Click the logout tab or “Logout” button in the home tab to logout from the app.



By default you need to sign in only when you start the application for the first time. If you are worried that somebody might get access to your files if your Phone is lost or stolen, you can disable the automatic login on your Phone .In this case, you will have to login each time you start the application.

## 03. HOW TO


### 3.1 HOW TO DO THE EYE TEST

Click the “start” button to start the eye test.

Before start the testing ensure you away 10ft or 20ft from the device. You can select 10 ft testing and 20ft testing parts from the web application.


Keep the required distance between the user and device.

Then follow the following steps.

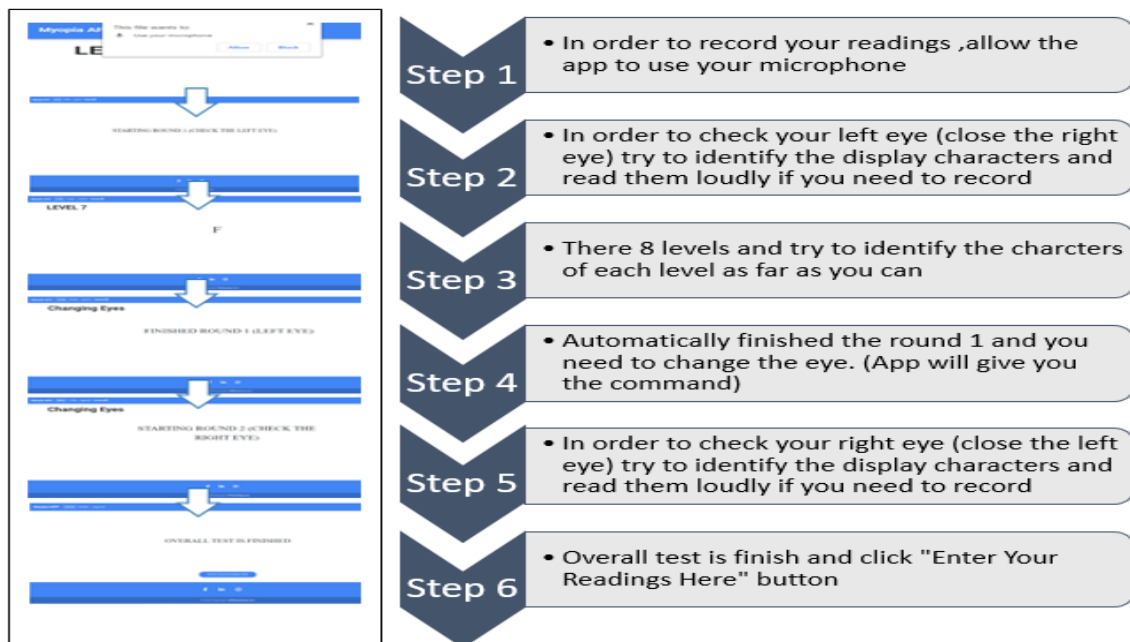


## 3.2 HOW TO CLARIFY THE RESULTS

Click “Enter Your Readings Here” button to create the results report.

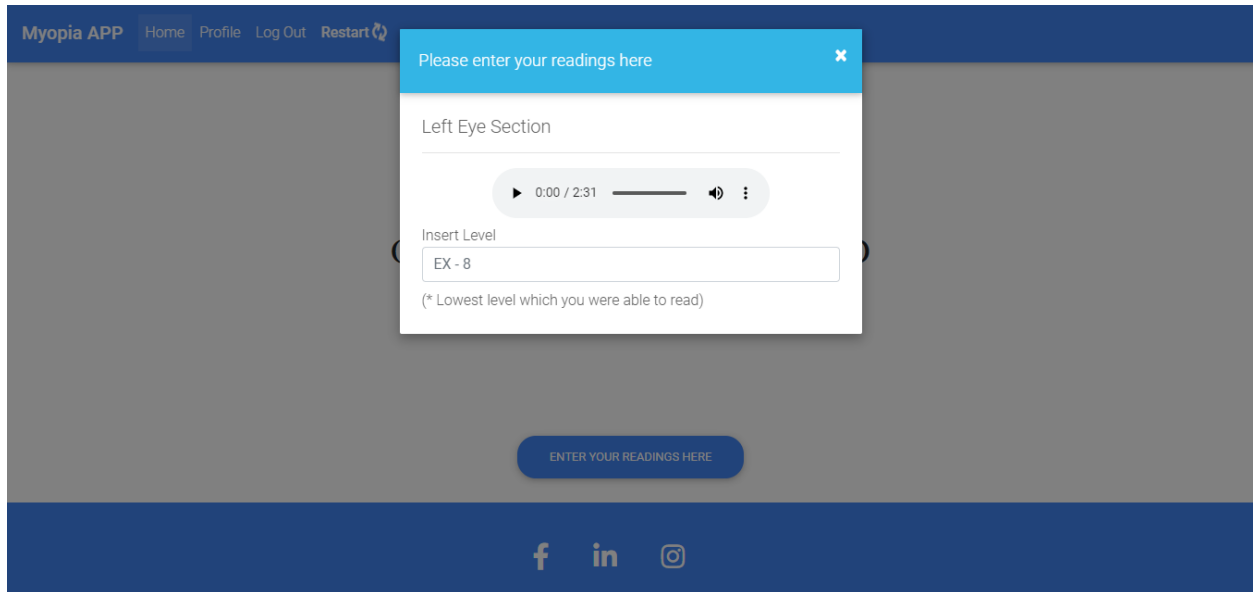

Here you can listen to your recordings and it allows you to find the last level and its characters which were identified by you.

Also there is an option to download your voice recordings as your preference.

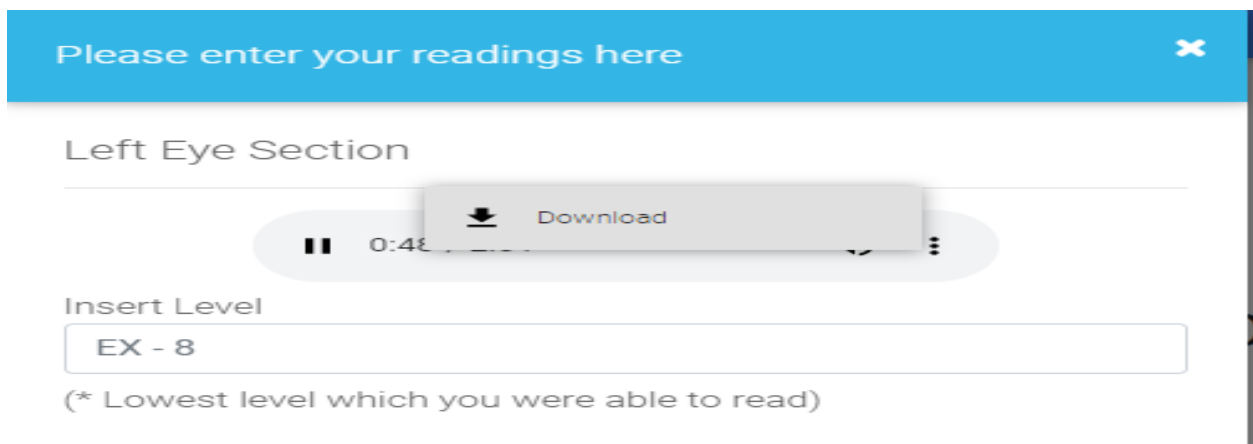

Enter the level and characters which were identified by your left eye at last.

Click the “Left accuracy” button to view your left eye accuracy.

Do the same in order to view your right eye accuracy.

Please enter your readings here

Left Eye Section

▶ 2:31 / 2:31

Insert Level

7

(\* Lowest level which you were able to read)

Insert characters

AKDC

(\* Don't use commas to separate characters)

LEFT ACCURACY

Left EYE Accuracy

(20/25)-7

Right Eye Section

▶ 0:00 / 5:07

Insert Level

8

(\* Lowest level which you were able to read)

Insert characters

BNHG

(\* Don't use commas to separate characters)

RIGHT ACCURACY

Right EYE Accuracy



(20/20)-8

CONFIRM

Click the “Confirm” button to confirm your data and then you may look an interface which allows you to edit your data or allows you to go the next page.



Press 'NEXT' and see your overarll Testing Results

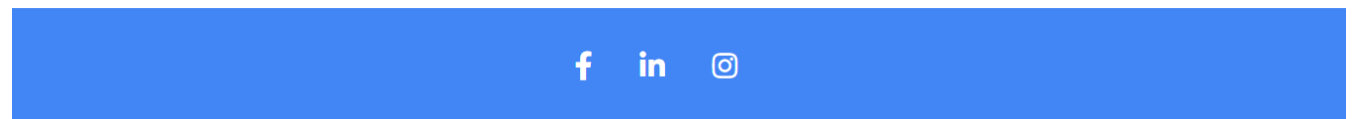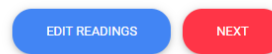

Click “Next” button to view your final results report and then you have an option to submit the data.

FINAL RESULT

LEFT EYE READINGS

Lowest Level you were able to recognaized

7

Orginal characters appeared

O,Z,T,P,L,E,Z

Your readings

A,K,Q,C

Accuracy

(20/25)-7

RIGHT EYE READINGS

Lowest Level you were able to recognaized

8

Orginal characters appeared

E,Z,Z,P,C,C,T,T

Your readings

B,N,H,G

Accuracy

(20/20)-8

SUBMIT DATA

## Do you know what does this mean?

Accuracy

(20/25)-7

A visual acuity of 20/25 is frequently described as meaning that a person can see detail from 20ft away the same as a person with "normal" eyesight would see from 25 ft. If a person has a visual acuity of (20/25)-7, you are said to see detail from 20 ft away the same as a person with "normal" eyesight would see it from 25 ft away. But in this current status you had missed 7 characters in the snellen chart.